


PPES-Series Pharmaceutical Grade Polyethersulfone

PPES-Series High Purity Pharmaceutical Grade
Polyethersulfone Filter Cartridges are ideal for sterile
filtration and clarification of pharmaceutical and
biological solutions. Each PPES cartridge is integrity
tested during manufacturing and is supported by
a validation guide for regulatory compliance. Low
protein binding and the broad chemical compatibility
characteristics of the polyethersulfone membrane,
along with exceptional flow rate vs pressure drop,
makes the PPES-Series the ideal choice for a
variety of valuable and/or critical pharmaceutical
solutions.

PPES cartridges are fully validated as sterilizing grade filters in accordance with HIMA and ASTM F838-05 guidelines. For the 0.2 micron series elements, validation studies demonstrate sterile effluent is achieved with challenge concentrations in excess of 10⁷, *Brevundimonas diminuta* per cm² of filter area. Additionally, validation studies of 0.1 micron series elements demonstrate 10⁷ retention of Mycoplasma (*Acholeplasma laidlawii*) per cm² of filter area. Designed to tolerate repeated hot water sanitization and *in-situ* steam sterilization cycles for maximum service life. Manufactured in a clean-room environment to maintain high standards of purity and cleanliness.



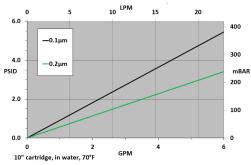



Typical Applications

- Vaccines
- Large Volume
 Parenteral (LVP's)
- Water for Injection (WFI)
- Diagnostics
- Ophthalmics
- Cell and Tissue Culture Media
- Protein Solutions
- Serum and Blood Products

Construction Materials

Membrane	Polyethersulfone
Support Media	Polypropylene
End Caps	Polypropylene
Center Core	Polypropylene
Outer Support Cage	Polypropylene
O-Rings/Gaskets	Buna, EPDM, Silicone,
Teflon® Encapsula	ated Viton®, Viton®, Teflon®
	Encapsulated Silicone

Note: O-ring adapters include integral reinforcement that will not deform with repeated steam sterilization or hot water sanitation cycles.

Toxicity

All polypropylene components meet the specifications for biological safety per USP Class VI – 121°C for plastics.

Sterilization

Hot Water	85°- 95°C, 30 min.,	max. ∆P 7 ps
In-Line Steaming] .13	34°C, 30 min.
	max. ∆P 7 p	si; 100 cycle

Dimensions

Length:

10 to 40 inches (25.4 to 101.6 cm) nominal

Outside Diameter:

2.78 inches (7.06 cm) nominal

Operating Conditions

Change Out ΔP (recommer	nded 35 PSID
Temperature (max)	176°F (80°C)
Differential Pressure (max)	72 PSID
	(5.0 bar) at 68°F (20°C)

Food Safety Compliance

Materials of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulations, 21CFR. Materials used to produce filter media and hardware are deemed safe for use in contact with foodstuffs in accordance with EU Directives 2002/72/EC, 1935/2004, and/or 10/2011.

Ordering Information

PPES	Rating (µ)	Α	Length	С	End Cap Style	O-Rings/Gaskets
	0.1		10" (25.4 cm)		2 = DOE Flat Gasket	B = Buna-N
	0.2		20" (50.8 cm)		3 = 222 w/ Fin	E = EPDM
			30" (76.2 cm)		4 = 222 w/ Flat Cap	S = Silicone
			40" (101.6 cm)		6 = 226 w/ Flat Cap	T = Teflon® Encapsulated Viton®
					7 = 226 w/ Fin	V = Viton®
					28 = 222 3-tabs w/ Fin	Z = Teflon® Encapsulated Silicone

DISCLAIMER: Filtration data presented is representative of performance observed in controlled laboratory testing. It is not given as a warranty, specification or statement of fitness for use. Specific performance can vary widely depending on contaminant type, fluid properties, flow rates and environmental conditions. It is recommended that users conduct thorough qualification testing to assure the product functions as required. For additional technical support, a product Validation Guide is available upon request.

DS PPES 190918

